

Daily Cleaning

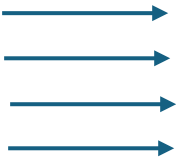
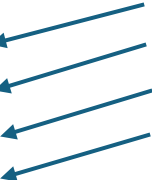
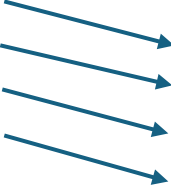
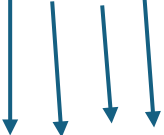
Before work:
Put on lab shoes before entering the lab
Put on lab coat and gloves
For bench tops: Clean workspace with 1:10 bleach dilution (DO NOT use bleach directly on metal or glass surfaces) Wait two minutes Rinse workspace with distilled water Clean workspace with 70% Ethanol
For LAF cabinets and fume hoods: Turn on air and UV sterilization for 5-10 minutes Spray low-lint cloth with 1:10 bleach dilution and wipe work surface back to front in straight lines Wait two minutes Rinse surface with DI water and wipe with low-lint cloth Spray low-lint cloth with 70% EtOH and wipe work surface
Clean pipettes and racks you will use with 70% Ethanol
Put clean waste bag on bench
After work:
Clean pipettes and racks with ethanol
Repeat cleaning for bench tops or cabinets
Empty bench waste bags into larger bins (plastic if it contains only tubes and tips or residual waste if it is mixed)
Remove lab coat and gloves
Bring full waste bags and cardboard boxes to outdoor trash and recycling
Bring solid chemical waste or full liquid chemical waste bottles to waste shed
Remove lab shoes outside the lab

Weekly Cleaning

Plan to spend 1 hour.

Task		Initial and date
Tidy up the room	☐	
Mix new 1:10 bleach dilution	☐	
Refill 70% EtOH bottles and DI water	☐	
Refill paper towels, gloves, tips	☐	
Wipe down the following with 1:10 bleach and 70% EtOH (DO NOT use bleach directly on metal or glass)	☐	
Countertops	☐	
Wipe down the following with 70% EtOH	☐	
Door, freezer, oven and cabinet handles	☐	
Pipettes	☐	
Chairs	☐	
Tops of trash bins and dusty surfaces	☐	
Sinks	☐	
Discard any sitting buffer in gel electrophoresis boxes	☐	
Weekly clean LAF cabinets and fume hoods	☐	
Wipe down machines with dry lint-free towel	☐	
Clean up storage room and empty plastic waste	☐	
Bring full waste bags and cardboard boxes to outdoor trash and recycling	☐	

Deep Clean LAF cabinets

- Always use low-lint towels and wipe in straight lines!
- Do not wipe the filter with any chemicals!
- Wet low-lint towel with 1:10 Bleach dilution until damp. Replace as needed.
 - Wipe the back left to right starting at the top and working your way down
 - 
 - Wipe the sides from back to front starting at the top and working your way down
 - 
 - 
 - Clean the work surface back to front, starting at the left and moving right
 - 
- Rinse with water
- Repeat with 70% Ethanol

Date last cleaned	Initials

Guidelines for moving objects from High-DNA/Post-PCR to Low-DNA/Pre-PCR space

1. In the post-PCR lab, wash item with 1:10 bleach dilution, rinse with water and wipe with 70% Ethanol solution
2. Bring to Museomics office space 05.2.028 and let sit overnight
3. The following day, bring the item to pre-PCR washroom and wash with bleach, water and ethanol
4. Your item is now ready to use in pre-PCR space!

Note: Autoclaved items may be brought into pre-PCR space if you wish to skip the overnight wait

Cleaning sign up

- Each lab user is responsible for cleaning once a month on Tuesday. Cleaning should take 1 hour. If you cannot attend your time, you need to switch times with someone else or find a replacement.

[illegible]